



賽馬會復康用品租賃中心

**Jockey Club Rehabilitation Products Rental Centre**

復康用品目錄

Rehabilitation Products Catalogue

**賽馬會復康用品租賃中心**為社區上有需要人士，提供完備及高質素的復康用品租賃服務。由用品選配、租賃、送遞、回收、全面清潔消毒及重新包裝等一應俱全。

### 服務對象

長者及社區上有需要人士

### 服務特色

- 專業評估/建議：設有義肢矯形師或職業治療師提供評估及諮詢服務，為顧客建議及選配合適復康用品
- 租金靈活彈性：租賃期及租金每星期(7天)計算，切合顧客不同需要
- 衛生安全管理：所有回收的復康用品均作全面清潔消毒、檢查維修及重新包裝，確保顧客能安心使用

### 租而不買的四大理由

- 支援短暫需要：切合短暫需要，如：復康用品故障、維修、購置中，幫助使用者繼續在社區生活
- 迎合復康階段：使用者可按身體狀況及需要轉變，於不同復康階段租用不同用品，促進復康過程
- 先租用後購買：復康用品種類及功能繁多，租借可透過親身試用，幫助使用者選取最合適的復康用品
- 推動環境保護：減少閒置、浪費與丟棄復康用品，符合「循環經濟、減少浪費、友善環境」的環保精神

### 服務收費

- 免費評估及諮詢
- 復康用品租金以每星期(7天)計算，視乎不同產品而定
- 租用人需繳付按金及運輸費(如適用)

## Rental Product Information 租賃產品資料

<b>Wheelchair 輪椅</b> .....	<b>6</b>
NISSIN L7 .....	6
NISSIN CBW .....	7
KARMA KM-5000* .....	8
CONVAID CX-12*, CX-14*, CX-16* .....	9
<b>Wheelchair with Tilting 可後傾座位輪椅</b> .....	<b>10</b>
NISSIN NAH-F5 (16" & 18") .....	10
KARMA VIP515 (16" & 18")* .....	11
CONVAID RD-12* & RD-14* .....	12
MARLEY RMS SIZE 1* & SIZE 2* .....	13
<b>Power Wheelchair 電動輪椅</b> .....	<b>14</b>
KD SUKD-00003 .....	14
MASAR MA-73 .....	15
OTTOBOCK JUVO B3 .....	16
<b>Power Wheelchair with Tilting 可後傾座位電動輪椅</b> .....	<b>17</b>
卓越 SUPREME VE7900 .....	17
<b>Wheelchair with Adaptive Seating 配置塑形座位輪椅</b> .....	<b>18</b>
KI MOBILITY CATALYST 4C* .....	18
SMART12*, SMART14* .....	19

<b>Wheelchair with Tilting Adaptive Seating 配置可後傾塑形座位輪椅</b> .....	<b>20</b>
KI MOBILITY FOCUS CR* .....	20
KI MOBILITY FLIP 14* .....	21
KI MOBILITY LIBERTY* .....	22
<b>Commode Chair for Adult with Tilting 可後傾座位浴椅</b> .....	<b>23</b>
SHOWERBUDDY SB3T* .....	23
AQUATEC OCEAN VIP COMMmode .....	24
<b>Commode Chair for Kid with Tilting 可後傾座位小童浴椅</b> .....	<b>25</b>
RIFTON HTS-Z120* .....	25
TUMBLE FORMS STARFISH BATH CHAIR SIZE3* .....	26
SHOWERBUDDY SB3T (FOR KID)* .....	27
<b>Hospital Bed (Four Folds) 電動復康床(四摺)</b> .....	<b>28</b>
SIDHIL MAGIC SIDHIL-011* .....	28
CAPRI BED FOLDING FLOORLINE BED* .....	29
<b>Foam / Air Mattress for Pressure Relief 泡棉 / 充氣減壓床墊</b> .....	<b>30</b>
SUPREME TURNING BED CUBE 2700 .....	30
INVACARE SOFTFORM PREMIER.....	31
<b>Mini-Hoist 小型移位機</b> .....	<b>32</b>
FORTA MINI-HOIST* .....	32
<b>Appendix I 附錄一</b> .....	<b>33</b>
Wheelchair Dimension Diagram 輪椅尺寸量度圖 .....	33

<b>Appendix II 附錄二</b> .....	<b>34</b>
Wheelchair Parts Identification 輪椅部件名稱 .....	34
<b>Appendix III 附錄三</b> .....	<b>40</b>
What is Adaptive Seating? .....	40
甚麼是配置塑形座位?.....	40
<b>Appendix IV 附錄四</b> .....	<b>41</b>
What are the differences between Tilting and Recline functions? .....	41
座位後傾和座背後傾功能有什麼分別?.....	42
<b>服務流程</b> .....	<b>43</b>

# Wheelchair 輪椅

## NISSIN L7

Model 型號	L7
Foldable 摺合	Yes 可以*
Seat Width (inch) 座闊(寸)	16
Seat Depth (inch) 座深(寸)	15
Seat Height (inch) 座高(寸)	17
Back Height (inch) 背高(寸)	17
Overall Size 整體尺寸	
Depth (inch) 深(寸)	36
Height (inch) 高(寸)	35
Width (inch) 闊(寸)	22
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	35
Height (inch) 高(寸)	27**
Width (inch) 闊(寸)	11
Max. User Weight (Kg) 最高承重量 (公斤)	75
Net Weight (Kg) 重量(公斤)	7.8



\*Foldable backrest 可摺背

\*\* With folded backrest 配可摺合椅背

返回目錄 

# Wheelchair 輪椅

## NISSIN CBW

Model 型號	CBW
Foldable 摺合	Yes 可以*
Seat Width (inch) 座闊(寸)	18
Seat Depth (inch) 座深(寸)	16
Seat Height (inch) 座高(寸)	18.5
Back Height (inch) 背高(寸)	15
Overall Size 整體尺寸	
Depth (inch) 深(寸)	38
Height (inch) 高(寸)	34
Width (inch) 闊(寸)	24.5
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	34
Height (inch) 高(寸)	26**
Width (inch) 闊(寸)	12
Max. User Weight (Kg) 最高承重量 (公斤)	100
Net Weight (Kg) 重量(公斤)	14.9

\*Foldable backrest 可摺背

\*\*Wwith folded backrest 配可摺合椅背



返回目錄



# Wheelchair 輪椅

## KARMA KM-5000\*

Model 型號	KM-5000
Foldable 摺合	Yes 可以
Seat Width (inch) 座闊 (寸)	16
Seat Depth (inch) 座深 (寸)	17
Seat Height (inch) 座高 (寸)	20
Back Height (inch) 背高 (寸)	18
Overall Size 整體尺寸	
Depth (inch) 深(寸)	50
Height (inch) 高(寸)	37
Width (inch) 闊(寸)	23
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	36
Height (inch) 高(寸)	37
Width (inch) 闊(寸)	13
Max. User Weight (Kg) 最高承重量 (公斤)	115
Net Weight (Kg) 重量 (公斤)	20.6
Recline 座背後傾角度	5°- 75°



\*With headrest detached 配可拆除頭枕

[返回目錄](#)





## Wheelchair 輪椅

### CONVAID CX-12\*, CX-14\*, CX-16\*

Model 型號	CX-12	CX-14	CX-16
Foldable 摺合	Yes 可以	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	12	14	16
Seat Depth (inch) 座深(寸)	8 – 13	10 – 15	11 – 16
Seat Height (inch) 座高(寸)	25	21	23
Back Height (inch) 背高(寸)	22	25	27
Overall Size 整體尺寸			
Depth (inch) 深(寸)	44	49	52
Height (inch) 高(寸)	38	37	40
Width (inch) 闊(寸)	23.5	25.5	27.5
Folded Size 摺合後尺寸			
Depth (inch) 深(寸)	48	47	53
Height (inch) 高(寸)	17	17	17
Width (inch) 闊(寸)	15	16	16
Max. User Weight (Kg) 最高承重量 (公斤)	30	45.5	77
Net Weight (Kg) 重量(公斤)	12	13	14
Fixed Seat Tilt 固定座位後傾度	30°	30°	30°



\*可配備太陽擋、加高頭墊、頭側翼、肩帶、可調校軟腰墊、分腿帶、小腿擋帶、腳帶、及/或座椅下儲物袋

\*Parts / accessories available : Canopy, Headrest extension, Padded Headwings, H-Harness, Double Flap Soft Adjustable Lateral Support with Scoliosis Strap, Medial Thigh Support Abductor, Calf Panel, Foot Positioners, and/or Under Seat Storage Basket

返回目錄



## Wheelchair with Tilting 可後傾座位輪椅

### NISSIN NAH-F5 (16" & 18")

Model 型號	NAH-F5 16"	NAH-F5 18"
Foldable 摺合	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	16	18
Seat Depth (inch) 座深(寸)	16	16
Seat Height (inch) 座高(寸)	16	16
Back Height (inch) 背高(寸)	19	19
Overall Size 整體尺寸		
Depth (inch) 深(寸)	43	43
Height (inch) 高(寸)	40 – 50	40 – 50
Width (inch) 闊(寸)	22	24
Folded Size 摺合後尺寸		
Depth (inch) 深(寸)	41	41
Height (inch) 高(寸)	26	26
Width (inch) 闊(寸)	13	13
Max. User Weight (Kg) 最高承重量 (公斤)	100	100
Net Weight (Kg) 重量(公斤)	20	20
Seat Tilt 座位後傾角度	-4° – 30°	-4° – 30°



返回目錄




## Wheelchair with Tilting 可後傾座位輪椅

### KARMA VIP515 (16" & 18")\*

Model 型號	VIP-515 16"	VIP-515 18"
Foldable 摺合	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	16	18
Seat Depth (inch) 座深(寸)	17	17
Seat Height (inch) 座高(寸)	19	19
Back Height (inch) 背高(寸)	18	18
Overall Size 整體尺寸		
Depth (inch) 深(寸)	43	43
Height (inch) 高(寸)	39	39
Width (inch) 闊(寸)	24.5	26
Folded Size 摺合後尺寸		
Depth (inch) 深(寸)	29	29
Height (inch) 高(寸)	30*	30*
Width (inch) 闊(寸)	13	13
Max. User Weight (Kg) 最高承重量 (公斤)	115	115
Net Weight (Kg) 重量(公斤)	18.9	19.8
Seat Tilt 座位後傾角度	0°- 35°	0°- 35°



\*With headrest detached 配可拆除頭枕

返回目錄 

# Wheelchair with Tilting 可後傾座位輪椅


## CONVAID RD-12\* & RD-14\*

Model 型號	RD-12	RD-14
Foldable 摺合	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	12	14
Seat Depth (inch) 座深(寸)	9 – 14	11 – 16
Seat Height (inch) 座高(寸)	21	22
Back Height (inch) 背高(寸)	22	24
Overall Size 整體尺寸		
Depth (inch) 深(寸)	35	35
Height (inch) 高(寸)	43	44
Width (inch) 闊(寸)	23	25
Folded Size 摺合後尺寸		
Depth (inch) 深(寸)	34	34
Height (inch) 高(寸)	17	18
Width (inch) 闊(寸)	23	25
Max. User Weight (Kg) 最高承重量 (公斤)	30	45.5
Net Weight (Kg) 重量(公斤)	15	17
Seat Tilt 座位後傾角度	5° – 45°	5° – 45°
Recline 座背後傾角度	90°/ 110°	90°/ 110°



\*可配備太陽擋、加高頭墊、頭側翼、肩帶、側墊、檯面板、分腿帶、小腿擋帶、腳帶、及/或座椅下儲物袋

\*Parts / accessories available : Canopy, Headrest extension, Padded Headwings, H-Harness, Lateral Contoured Trunk Support, Tray, Medial Thigh Support Abductor, Calf Panel, Foot Positioners, and/or Under Seat Storage Basket

返回目錄 

## Wheelchair with Tilting 可後傾座位輪椅


### MARLEY RMS SIZE 1\* & SIZE 2\*

Model 型號	RMS Size 1	RMS Size 2
Foldable 摺合	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	8 – 11	10 – 13
Seat Depth (inch) 座深(寸)	9 – 11	11 – 14
Seat Height (inch) 座高(寸)	19	20.5
Back Height (inch) 背高(寸)	20	22
Overall Size 整體尺寸		
Depth (inch) 深(寸)	31	32
Height (inch) 高(寸)	42	46
Width (inch) 闊(寸)	23	25
Folded Size 摺合後尺寸		
Depth (inch) 深(寸)	31.5	33
Height (inch) 高(寸)	20	20
Width (inch) 闊(寸)	23	25
Max. User Weight (Kg) 最高承重量 (公斤)	32	32
Net Weight (Kg) 重量(公斤)	14	16
Seat Tilt 座位後傾角度	15°–50°	15°–50°



\*可配備太陽擋、加高頭墊、頭側翼、肩帶、側墊、檯面板、腳帶、及/或座椅下儲物袋

\*Parts / accessories available : Canopy, Headrest extension, Padded Headwings, H-Harness, Lateral Trunk Support, Tray, Foot Positioners, and/or Under Seat Storage Basket


返回目錄 

# Power Wheelchair 電動輪椅

## KD SUKD-00003

Model 型號	SUKD-00003
Foldable 摺合	Yes 可以
Seat Width (inch) 座闊(寸)	15
Seat Depth (inch) 座深(寸)	17
Seat Height (inch) 座高(寸)	18.5
Back Height (inch) 背高(寸)	14
Overall Size 整體尺寸	
Depth (inch) 深(寸)	39
Height (inch) 高(寸)	34
Width (inch) 闊(寸)	26
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	13
Height (inch) 高(寸)	30
Width (inch) 闊(寸)	26
Max. User Weight (Kg) 最高承重量 (公斤)	115
Net Weight (Kg) 重量(公斤)	23




返回目錄 

# Power Wheelchair 電動輪椅

## MASAR MA-73

Model 型號	Ma-73
Foldable 摺合	Yes 可以
Seat Width (inch) 座闊(寸)	18
Seat Depth (inch) 座深(寸)	17
Seat Height (inch) 座高(寸)	19
Back Height (inch) 背高(寸)	15
Overall Size 整體尺寸	
Depth (inch)	37
Height (inch)	35
Width (inch)	26
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	37
Height (inch) 高(寸)	35
Width (inch) 闊(寸)	20
Max. User Weight (Kg) 最高承重量 (公斤)	125
Net Weight (Kg) 重量(公斤)	22



返回目錄 

# Power Wheelchair 電動輪椅

## OTTOBOCK JUVO B3

Model 型號	Juvo B3 16"	Juvo B3 18"
Foldable 摺合	No 不可以	No 不可以
Seat Width (inch) 座闊(寸)	16	18
Seat Depth (inch) 座深(寸)	16	18
Seat Height (inch) 座高(寸)	19	20
Back Height (inch) 背高(寸)	22	22
Overall Size 整體尺寸		
Depth (inch) 深(寸)	42	42
Height (inch) 高(寸)	42	42
Width (inch) 闊(寸)	23.5	23.5
Max. User Weight (Kg) 最高承重量 (公斤)	113	113
Net Weight (Kg) 重量(公斤)	93	93
Recline 座背後傾角度	90°/ 100°/110° /120°	90°/ 100°/110° /120°



返回目錄 




## Power Wheelchair with Tilting 可後傾座位電動輪椅

### 卓越 SUPREME VE7900

Model 型號	VE7900
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	16
Seat Depth (inch) 座深(寸)	18
Seat Height (inch) 座高(寸)	21
Back Height (inch) 背高(寸)	21
Overall Size 整體尺寸	
Depth (inch) 深(寸)	39
Height (inch) 高(寸)	49
Width (inch) 闊(寸)	25
Max. User Weight (Kg) 最高承重量 (公斤)	136
Net Weight (Kg) 重量(公斤)	115
Seat Tilt 座位後傾角度	3° - 21°



返回目錄 

# Wheelchair with Adaptive Seating 配置塑形座位輪椅

## KI MOBILITY CATALYST 4C\*

Model 型號	CATALYST 4C 14"	CATALYST 4C 16"	CATALYST 4C 18"
Foldable 摺合	Yes 可以	Yes 可以	Yes 可以
Seat Width (inch) 座闊(寸)	14	16	18
Seat Depth (inch) 座深(寸)	14	16/18	16
Seat Height (inch) 座高(寸)	16.5	16	16
Back Height (inch) 背高(寸)	18.5	17.5	17.5
Overall Size 整體尺寸			
Depth (inch) 深(寸)	35	39	39
Height (inch) 高(寸)	37	36	36
Width (inch) 闊(寸)	23	24	26
Folded Size 摺合後尺寸			
Depth (inch) 深(寸)	35	39	43
Height (inch) 高(寸)	37	36	36
Width (inch) 闊(寸)	14	14	14
Max. User Weight (Kg) 最高承重量 (公斤)	113	113	113
Net Weight (Kg) 重量(公斤)	11	11	11



\*可配備頭托、泡棉背墊、側墊、檯面板、高回彈泡棉座墊、點滴架、防翻桿、小腿擋帶及/或可升降腳踏

\*Parts / accessories available: Headrest, Backrest with Foam Paddings, Lateral Pads, Tray, Viscoelastic Seat Cushion, IV Pole, Anti-Tippers, Calf Straps and/or Elevating Leg Rest

返回目錄 

# Wheelchair with Adaptive Seating 配置塑形座位輪椅

## SMART12\*, SMART14\*

Model 型號	Smart12	Smart14
Foldable 摺合	No 不可以	No 不可以
Seat Width (inch) 座闊(寸)	12	14
Seat Depth (inch) 座深(寸)	10	12
Seat Height (inch) 座高(寸)	23	23
Back Height (inch) 背高(寸)	12 – 18	12 – 18
Overall Size 整體尺寸		
Depth (inch) 深(寸)	35	35
Height (inch) 高(寸)	47	47
Width (inch) 闊(寸)	22	22
Max. User Weight (Kg) 最高承重量 (公斤)	60	60
Net Weight (Kg) 重量(公斤)	24.6	24.6
Seat Tilt 座位後傾角度	0°– 40°	0°– 40°



\*可配備側墊、胸帶、分腿墊、防翻桿、及/或小腿擋帶

\* Parts / accessories available : Lateral Pads, Posture Support Vest, PU Abductor Pads, Anti-Tippers, and/or Calf Straps

返回目錄



# Wheelchair with Tilting Adaptive Seating 配置可後傾塑形座位輪椅

## KI MOBILITY FOCUS CR\*

Model 型號	Focus CR
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	16/18
Seat Depth (inch) 座深(寸)	17
Seat Height (inch) 座高(寸)	16.5
Back Height (inch) 背高(寸)	22 – 26
Overall Size 整體尺寸	
Depth (inch) 深(寸)	51
Height (inch) 高(寸)	49
Width (inch) 闊(寸)	24/26
Max. User Weight (Kg) 最高承重量 (公斤)	136
Net Weight (Kg) 重量(公斤)	15
Seat Tilt 座位後傾角度	-5° – 50°
Recline 座背後傾角度	0° – 65°



\*可配備太陽擋、側墊、心口帶、檯面板、防翻桿、小腿擋帶及/或可升降腳踏

\*Parts / accessories available : Canopy, Lateral Pads, Postural Support Strapped Vest, Tray, Anti-Tippers, Calf Straps and/or Elevating Leg Rest

返回目錄



# Wheelchair with Tilting Adaptive Seating 配置可後傾塑形座位輪椅


## KI MOBILITY FLIP 14\*

Model 型號	Flip 14
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	14
Seat Depth (inch) 座深(寸)	14
Seat Height (inch) 座高(寸)	16.5
Back Height (inch) 背高(寸)	20 – 22
Overall Size 整體尺寸	
Depth (inch) 深(寸)	40
Height (inch) 高(寸)	42
Width (inch) 闊(寸)	20
Max. User Weight (Kg) 最高承重量 (公斤)	75
Net Weight (Kg) 重量(公斤)	9.5
Seat Tilt 座位後傾角度	0°– 45°



\*可配備太陽擋、側墊、肩帶、防翻桿、及/或小腿擋帶

\*Parts / accessories available : Canopy, Lateral Pads, Shoulder Harness, Anti-Tippers, and/or Calf Straps

返回目錄 

# Wheelchair with Tilting Adaptive Seating 配置可後傾塑形座位輪椅


## KI MOBILITY LIBERTY\*

Model 型號	Liberty FT
Foldable 摺合	Yes 可以
Seat Width (inch) 座闊(寸)	18
Seat Depth (inch) 座深(寸)	18
Seat Height (inch) 座高(寸)	16
Back Height (inch) 背高(寸)	18
Overall Size 整體尺寸	
Depth (inch) 深(寸)	38
Height (inch) 高(寸)	40
Width (inch) 闊(寸)	26
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	43
Height (inch) 高(寸)	40
Width (inch) 闊(寸)	14
Max. User Weight (Kg) 最高承重量 (公斤)	113
Net Weight (Kg) 重量(公斤)	11.7
Seat Tilt 座位後傾角度	0°– 20°



\*可配備頭托、側墊、心口帶、防翻桿、及/或小腿擋帶

\*Parts / accessories available : Headrest, Lateral Pads, Postural Support Vest , Anti-Tippers, and/or Calf Straps

返回目錄 

## Commode Chair for Adult with Tilting 可後傾座位浴椅

### SHOWERBUDDY SB3T\*

Model 型號	sb3t
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	18
Seat Depth (inch) 座深(寸)	17
Seat Height (inch) 座高(寸)	20 – 24
Back Height (inch) 背高(寸)	17
Overall Size 整體尺寸	
Depth (inch) 深(寸)	43.3
Height (inch) 高(寸)	43 – 55
Width (inch) 闊(寸)	22.5
Max. User Weight (Kg) 最高承重量 (公斤)	159
Net Weight (Kg) 重量(公斤)	22
Seat Tilt 座位後傾角度	0° – 30°



\*可配備座墊蓋、及/或便兜

\*Parts / accessories available : Seat Cushion Overlay, and/or Pan

返回目錄




# Commode Chair for Adult with Tilting 可後傾座位浴椅

## AQUATEC OCEAN VIP COMMODE

Model 型號	AQUATEC OCEAN VIP COMMODE
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	20
Seat Depth (inch) 座深(寸)	18
Seat Height (inch) 座高(寸)	20 – 24
Back Height (inch) 背高(寸)	21
Overall Size 整體尺寸	
Depth (inch) 深(寸)	42
Height (inch) 高(寸)	48
Width (inch) 闊(寸)	23.5
Max. User Weight (Kg) 最高承重量 (公斤)	150
Net Weight (Kg) 重量(公斤)	23.5
Seat Tilt 座位後傾角度	0° – 35°



返回目錄 



## Commode Chair for Kid with Tilting 可後傾座位小童浴椅


### RIFTON HTS-Z120\*

Model 型號	HTS-Z120
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	10/ 12/14
Seat Depth (inch) 座深(寸)	11.5 – 14.5
Seat Height (inch) 座高(寸)	15.5 – 21.5
Back Height (inch) 背高(寸)	19.5
Overall Size 整體尺寸	
Depth (inch) 深(寸)	31
Height (inch) 高(寸)	46
Width (inch) 闊(寸)	21
Max. User Weight (Kg) 最高承重量 (公斤)	68
Net Weight (Kg) 重量(公斤)	12.5
Seat Tilt 座位後傾角度	-5° – 25°
Recline 座背後傾角度	84°/ 96°/108°



\*可配備推手、心口帶、側墊、檯面板、側擋、導流器、便兜、小腿擋帶、及/或腳帶

\*Parts / accessories available : Push Handles, Butterfly-Harness, Lateral Supports, Tray, Hip Guides, Deflectors, Pan, Calf Rest, and/or Ankle Straps

返回目錄 

## Commode Chair for Kid with Tilting 可後傾座位小童浴椅


### TUMBLE FORMS STARFISH BATH CHAIR SIZE3\*

Model 型號	Size3
Foldable 摺合	Yes 可以
Seat Width (inch) 座闊(寸)	16
Seat Depth (inch) 座深(寸)	13
Seat Height (inch) 座高(寸)	5 – 17
Back Height (inch) 背高(寸)	18
Overall Size 整體尺寸	
Depth (inch) 深(寸)	29 – 52.5
Height (inch) 高(寸)	5 – 38
Width (inch) 闊(寸)	16
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	52.5
Height (inch) 高(寸)	5
Width (inch) 闊(寸)	16
Max. User Weight (Kg) 最高承重量 (公斤)	72
Net Weight (Kg) 重量(公斤)	5.2
Recline 座背後傾角度	90°–180°



\*可配備頭側翼、大腿帶、胸帶、側墊、及/或分腿帶

\*Parts / accessories available : Head Laterals, Hip Belt, Chest Belt, side pad and/or Abduction Strap

返回目錄 

# Commode Chair for Kid with Tilting 可後傾座位小童浴椅


## SHOWERBUDDY SB3T (FOR KID)\*

Model 型號	sb3t (available for kid)
Foldable 摺合	No 不可以
Seat Width (inch) 座闊(寸)	17
Seat Depth (inch) 座深(寸)	15
Seat Height (inch) 座高(寸)	21.5 – 25.5
Back Height (inch) 背高(寸)	4 – 19
Overall Size 整體尺寸	
Depth (inch) 深(寸)	43.3
Height (inch) 高(寸)	42.5 – 46.5
Width (inch) 闊(寸)	19
Max. User Weight (Kg) 最高承重量 (公斤)	159
Net Weight (Kg) 重量(公斤)	22
Seat Tilt 座位後傾角度	0° – 30°



\*可配備便兜

\*Parts / accessories available : Pan

返回目錄 

# Hospital Bed (Four Folds) 電動復康床(四摺)

## SIDHIL MAGIC SIDHIL-011\*

Model 型號	SIDHIL-011
Foldable 摺合	Yes 可以
Overall Size 整體尺寸	
Depth (inch) 深(寸)	84
Height (inch) 高(寸)	9.6 – 22
Width (inch) 闊(寸)	38
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	17
Height (inch) 高(寸)	45
Width (inch) 闊(寸)	38
Max. User Weight (Kg) 最高承重量 (公斤)	213
Net Weight (Kg) 重量(公斤)	90
Recline 背部升降高度	0° – 74°



\*租賃包括普通床墊一張

\*ONE regular mattress inclusive per one hospital bed rental

[返回目錄](#)



# Hospital Bed (Four Folds) 電動復康床(四摺)

## CAPRI BED FOLDING FLOORLINE BED\*

Model 型號	MC310
Foldable 摺合	Yes 可以
Overall Size 整體尺寸	
Depth (inch) 深(寸)	92
Height (inch) 高(寸)	42
Effective Life Range (inch) 有效升降延伸(寸)	5.1 – 28.5
Width (inch) 闊(寸)	36
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	18
Height (inch) 高(寸)	46
Width (inch) 闊(寸)	36
Max. User Weight (Kg) 最高承重量 (公斤)	225
Net Weight (Kg) 重量(公斤)	120
Recline 背部升降高度	0° – 65°



\*租賃包括普通床墊一張

\*ONE regular mattress inclusive per one hospital bed rental

[返回目錄](#)




# Foam / Air Mattress for Pressure Relief 泡棉 / 充氣減壓床墊

## SUPREME TURNING BED CUBE 2700

Model 型號	CUBE 2700
Foldable 摺合	Yes 可以
Overall Size 整體尺寸	
Depth (inch) 深(寸)	75
Height (inch) 高(寸)	4
Width (inch) 闊(寸)	31
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	17
Height (inch) 高(寸)	12
Width (inch) 闊(寸)	14
Max. User Weight (Kg) 最高承重量 (公斤)	120
Net Weight (Kg) 重量(公斤)	8




返回目錄 

## Foam / Air Mattress for Pressure Relief 泡棉 / 充氣減壓床墊

### INVACARE SOFTFORM PREMIER

Model 型號	SFPremier
Foldable 摺合	No 不可以
Overall Size 整體尺寸	
Depth (inch) 深(寸)	77.5
Height (inch) 高(寸)	6
Width (inch) 闊(寸)	34
Max. User Weight (Kg) 最高承重量 (公斤)	247.6
Net Weight (Kg) 重量(公斤)	13



返回目錄 

## Mini-Hoist 小型移位機

### FORTA MINI-HOIST\*

Model 型號	Forta Mini-Hoist
Foldable 摺合	Yes 可以
Overall Size 整體尺寸	
Depth (inch) 深(寸)	41
Height (inch) 高(寸)	70
Width (inch) 闊(寸)	35
Folded Size 摺合後尺寸	
Depth (inch) 深(寸)	17
Height (inch) 高(寸)	51
Width (inch) 闊(寸)	21
Max. User Weight (Kg) 最高承重量 (公斤)	130
Net Weight (Kg) 重量(公斤)	37



\*租賃包括吊兜一個

\*ONE Sling inclusive per one mini hoist rental

Sling Size (kg) 吊兜大小(承重量(公斤)) : M 中 (45 – 95) / S 小 (25 –50)











[返回目錄](#)





# Appendix I 附錄一

## Wheelchair Dimension Diagram 輪椅尺寸量度圖

尺寸量度	Dimension Measurement	圖例 Illustration	尺寸量度	Dimension Measurement	圖例 Illustration
座闊	Seat Width		背高	Back Height	
座深	Seat Depth		座高	Seat Height	
摺合後深	Folded Depth		整體總深	Overall Depth	
摺合後高	Folded Height		整體總高	Overall Height	
摺合後闊	Folded Width		整體總闊	Overall Width	

返回目錄 


## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

正方圖示一

- |                                 |                                  |
|---------------------------------|----------------------------------|
| 1 頭枕<br>Headrest                | 7 枱面板<br>Tray                    |
| 2 推手<br>Push Handle             | 8 安全帶<br>Pelvic Positioning Belt |
| 3 推手高低可調校掣<br>Adjustable Handle | 9 坐墊<br>Seat Cushion             |
| 4 背墊<br>Back Upholstery         | 10 腿跟擋帶<br>Heel Loop             |
| 5 側墊<br>Lateral Support Padding | 11 可升降腳踏<br>Elevating Legrest    |
| 6 心口帶<br>Vest                   | 12 吊點滴架<br>IV Pole               |



返回目錄 


## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

正方圖示二

- 1 太陽擋  
Canopy
- 2 頭枕  
Headrest
- 3 背墊  
Back Upholstery
- 4 側墊  
Lateral Support Padding
- 5 肩帶  
Shoulder Harness
- 6 安全帶  
Pelvic Positioning Belt
- 7 坐墊  
Seat Cushion
- 8 腿跟擋帶  
Heel Loop



返回目錄 

## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

後方圖示一

- 1 推手  
Push Handle
- 2 手柄剎  
Attendant Handbrake
- 3 調校座背後傾手掣  
Recline Handle
- 4 背柱  
Back Cane
- 5 防翻桿  
Anti-tipper
- 6 剎掣  
Wheel lock



返回目錄



## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

後方圖示二

- 1 推手  
Push Handle
- 2 調校座背後傾手掣  
Recline Handle
- 3 調校座位後傾腳掣  
Foot Tilt
- 4 防翻桿  
Anti-tipper
- 5 吊點滴架  
IV Pole



返回目錄



## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

側面圖示一

- 1 推手  
Push Handle
- 2 推手高低可調校掣  
Adjustable Handle
- 3 調校座背後傾手掣  
Recline Handle
- 4 後輪  
Rear Wheel
- 5 剎掣  
Wheel lock
- 6 前輪  
Caster
- 7 可昇降腳踏  
Elevating Legrest
- 8 可掀起扶手  
Flip Back Armrest
- 9 可昇降扶手  
Height Adjustable Armrest



返回目錄




## Appendix II 附錄二

### Wheelchair Parts Identification 輪椅部件名稱

側面圖示二

- 1 背柱  
Back Cane
- 2 手柄剎  
Attendant Handbrake
- 3 摺疊背柱手掣  
Foldable Backcane Handle
- 4 小腿擋帶  
Calf Strap
- 5 前轆  
Caster
- 6 後轆  
Rear Wheel

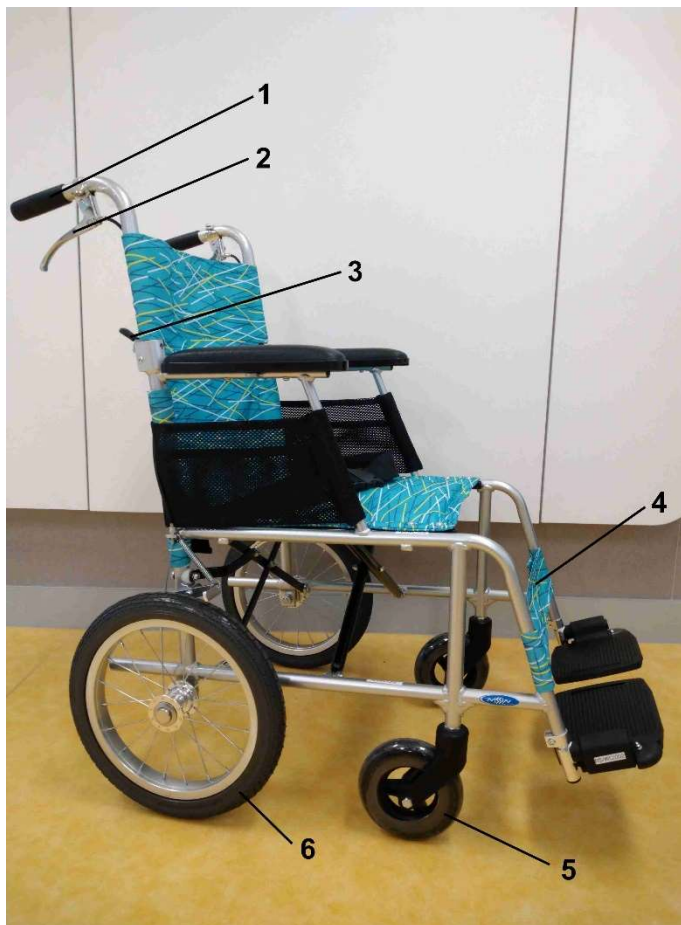

返回目錄



# What is Adaptive Seating?

## 甚麼是配置塑形座位?

- 配置塑形座位是透過義肢矯形師或治療師選擇不同設計配件的座椅，令使用者不但感到舒適、並會配合肢體變形、減低坐壓、改善其平衡和控制坐姿、以延長使用輪椅的時間，提升使用者於上學、工作、日常自理、社交、休閒玩樂或進行其他活動的能力。研究顯示，使用輪椅系統人士能從職能上感到滿足，相對下較樂於參與社交活動<sup>1</sup>。
- Adaptive Seating is tailor-made through the prescription of appropriate chair and/or its parts carefully chosen by a therapist or Prosthetist & Orthotist which are suitable for a user. It provides comfort, accommodates deformities, reduces pressure point, enhances user's sitting balance and posture, and extends sitting tolerance. Its goals are to improve wheelchair users' activities performance whether in school, at work, during daily activities, social participation and leisure, and other activities which are meaningful to them. Researchers observe that individuals who use adaptive seating devices are often gaining more occupational satisfaction, and thus increase their willingness to participate in social activities.
- 研究顯示，配置塑形座位能改善使用者的體位控制，從而提升使用者日常生活質，例如改善上肢的活動能力、抓握力、視覺掃描及追蹤能力；系統對提升自我進食的技巧尤其見效，原因是系統能維持使用者盤骨保持中線位置，讓脊髓有一個穩固的根基，以致上肢功能上引發的活動能力能有效發揮。總括而言，輪椅系統可促進軀幹的正常位置，令使用者能夠有餘力控制雙手及有助於消化食物<sup>2</sup>。
- Adaptive seating is believed to be able to improve postural control because it allows users to develop and enhance functional ability<sup>1,3</sup> , particularly mobility of the arms, grasping, visual scanning and tracking, and especially self-feeding. These devices are particularly useful in the development of self-feeding skills because they hold the pelvis in a neutral position which form a secure foundation for the rest of the spine, so as to allow the effective use of the upper limbs in functional activities. In conclusion, adaptive seating provides a proper trunk alignment and position which helps to improve the controlling of the hands and allows for proper digestion<sup>3</sup>.
- 配置塑形座位亦有效減低脊椎變型的風險。在某些情況下，它可以幫助矯正大腦麻痺小孩的脊骨側彎問題<sup>4</sup>。
- Adaptive seating is effective in reducing the risk of spinal deformity. In some cases, these devices can help in the correction of existing scoliosis in children with CP<sup>4</sup>.

[返回目錄](#)



<sup>1</sup> Rigby, P. J., Ryan, S. E., & Campbell, K. A. (2009). Effect of adaptive seating devices on the activity performance of children with cerebral palsy. *Archives of physical medicine and rehabilitation*, 90(8), 1389-1395.

<sup>2</sup> Hulme, J. B., Shaver, J., Acher, S., Mullette, L., & Eggert, C. (1987). Effects of adaptive seating devices on the eating and drinking of children with multiple handicaps. *American Journal of Occupational Therapy*, 41(2), 81-89.


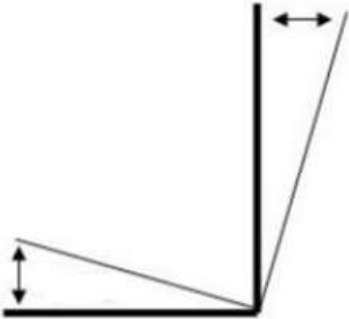

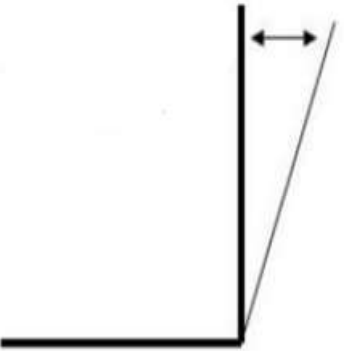
<sup>3</sup> Angsupaisal, M., Maathuis, C. G., & Hadders-Algra, M. (2015). Adaptive seating systems in children with severe cerebral palsy across International Classification of Functioning, Disability and Health for Children and Youth version domains: a systematic review. *Developmental Medicine & Child Neurology*, 57(10), 919-930.

<sup>4</sup> Holmes, K. J., Michael, S. M., Thorpe, S. L., & Solomonidis, S. E. (2003). Management of scoliosis with special seating for the non-ambulant spastic cerebral palsy population—a biomechanical study. *Clinical biomechanics*, 18(6), 480-487.



## Appendix IV 附錄四

### What are the differences between Tilting and Recline functions?

Tilting seat	Recline Seat
 	 
<ul style="list-style-type: none"> <li>• Tilting seat could tilt up the entire seat frame, keeping the seat-to-back angle unchanged as well as the positions of hip, knee, and ankle.</li> <li>• Designed for comfort, pressure relief and constant position; and it also reduces shearing force.</li> <li>• It relieves pressure beneath the pelvis when the seat is tilted backward for 15°; the risk of sliding out is further reduced when it is at 25°.</li> </ul>	<ul style="list-style-type: none"> <li>• Recline seats only tilt the back but not the seat as such the seat-to-back angle is changed as well as the lower extremity joint angles.</li> <li>• Seat to back angle at 120° would reduce pressure in ischial tuberosity but it does not decrease the pressure on sacrum. On the other hand, it might increase the chance of sliding out. Thus, an elevated footrest is required<sup>5</sup>.</li> </ul>


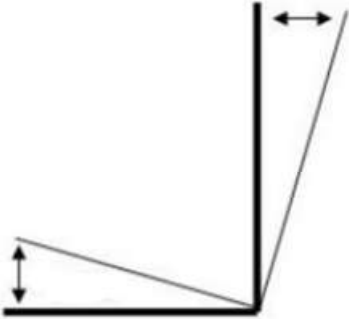

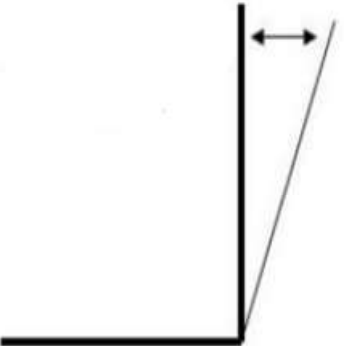
返回目錄



<sup>5</sup> Jan, Y. K., Crane, B. A., Liao, F., Woods, J. A., & Ennis, W. J. (2013). Comparison of muscle and skin perfusion over the ischial tuberosities in response to wheelchair tilt-in-space and recline angles in people with spinal cord injury. *Archives of physical medicine and rehabilitation*, 94(10), 1990-1996.

## Appendix IV 附錄四

### 座位後傾和座背後傾功能有什麼分別?


座位後傾	座背後傾
 	 
<ul style="list-style-type: none"><li>• 座位後傾功能指座位後傾時，座位和座背同時向後傾，此時座高增加，形同將整個輪椅座位向上翹起，而座椅與座背之間的角度不變，大腿、小腿和腳踝的位置亦保持不變。</li><li>• 座位後傾能減低盤骨下的壓力，減少剪力以改善坐姿。</li><li>• 座位後傾 15 度時能減少盤骨的壓力；座位後傾 25 度時，能減少向前滑出的危險。</li></ul>	<ul style="list-style-type: none"><li>• 座背後傾功能能增加座椅與座背之間的角度，但座高保持不變，下肢的角度則因應改變。</li><li>• 座背後傾 120 度時能減低坐骨的壓力，但無助減低骶骨的壓力，同時增加向前滑的風險，因此需要安裝可升降腳踏配合此功能<sup>5</sup>。</li></ul>

返回目錄



<sup>5</sup> Jan, Y. K., Crane, B. A., Liao, F., Woods, J. A., & Ennis, W. J. (2013). Comparison of muscle and skin perfusion over the ischial tuberosities in response to wheelchair tilt-in-space and recline angles in people with spinal cord injury. *Archives of physical medicine and rehabilitation*, 94(10), 1990-1996.

# 服務流程

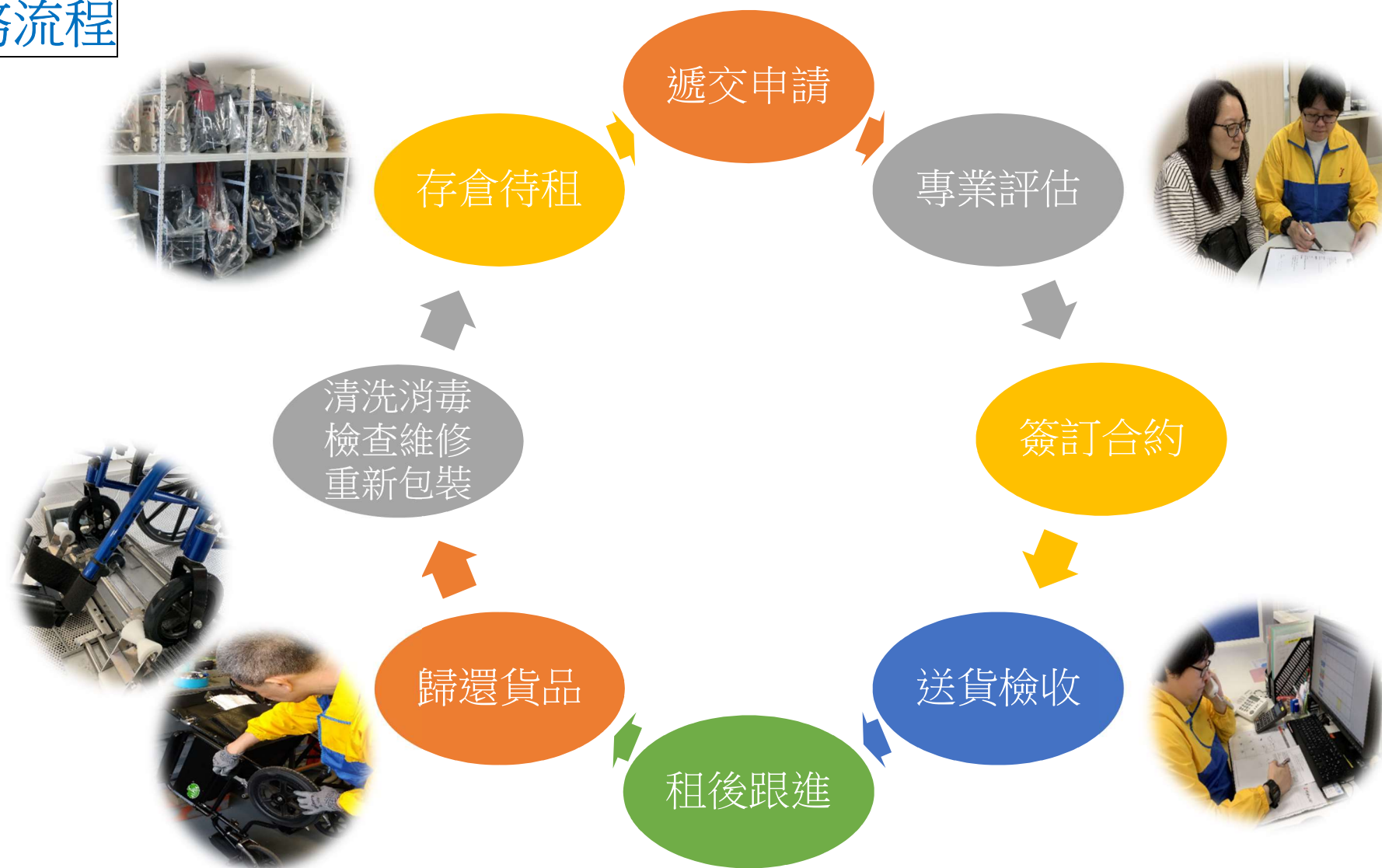

## 服務申請

- 向本中心職員索取《服務申請表》，填寫後交回本中心

[返回目錄](#)



地址：新界沙田火炭山尾街 5-13 號穗輝工廠大廈 11 樓 1101-1107 室

Address：Unit 1101-1107 11/F. Sui Fai Factory Estate 5-13 Shan Mei Street Fo Tan Shatin New Territories

電話 Tel：(852) 2736 8988

傳真 Fax：(852) 2617 3689

電郵 Email：[rental@sahk1963.org.hk](mailto:rental@sahk1963.org.hk)